





A perfect prosthesis begins with a  
perfect impression.....






## Edentulous Impression Trays

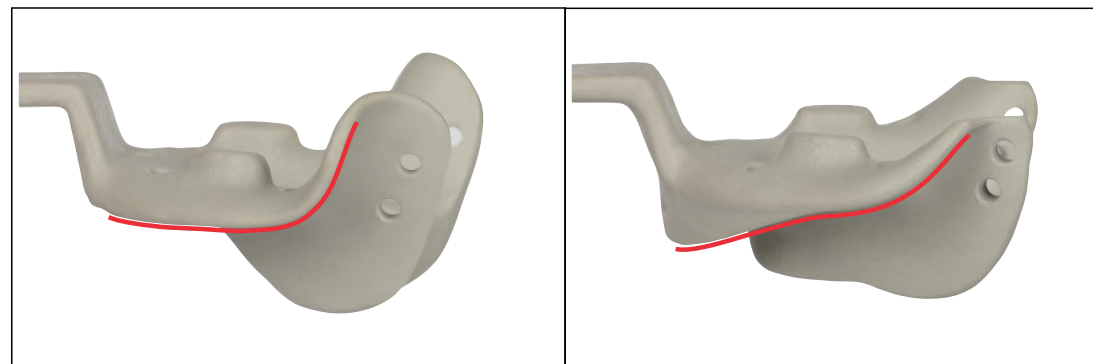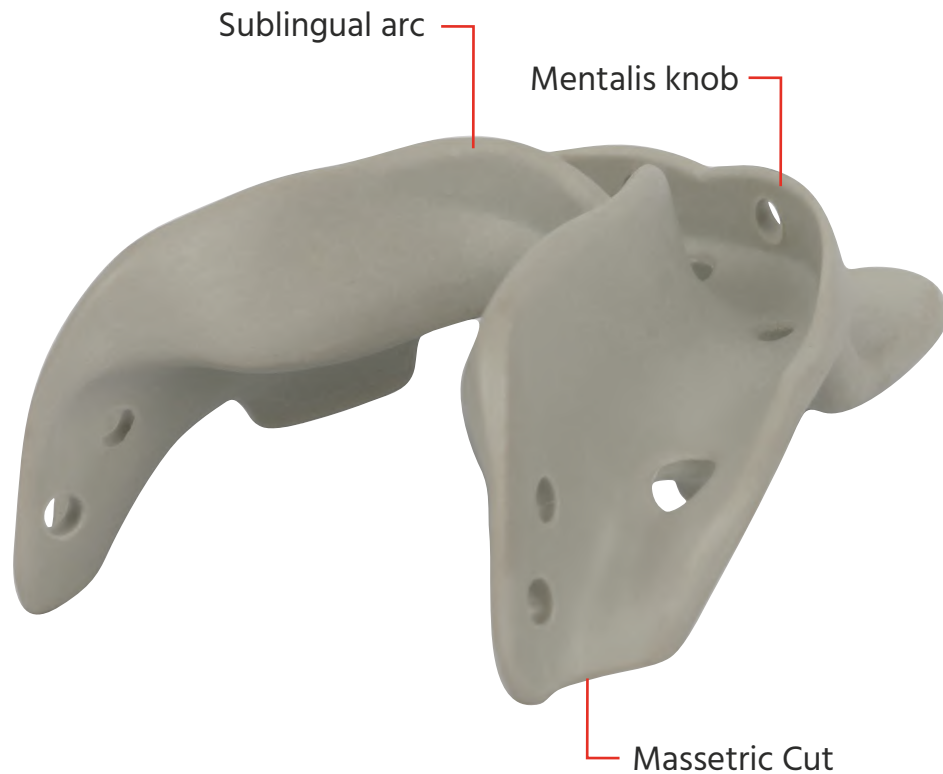


Tray Contour for Atrophied Ridge

Tray Contour for Average Ridge

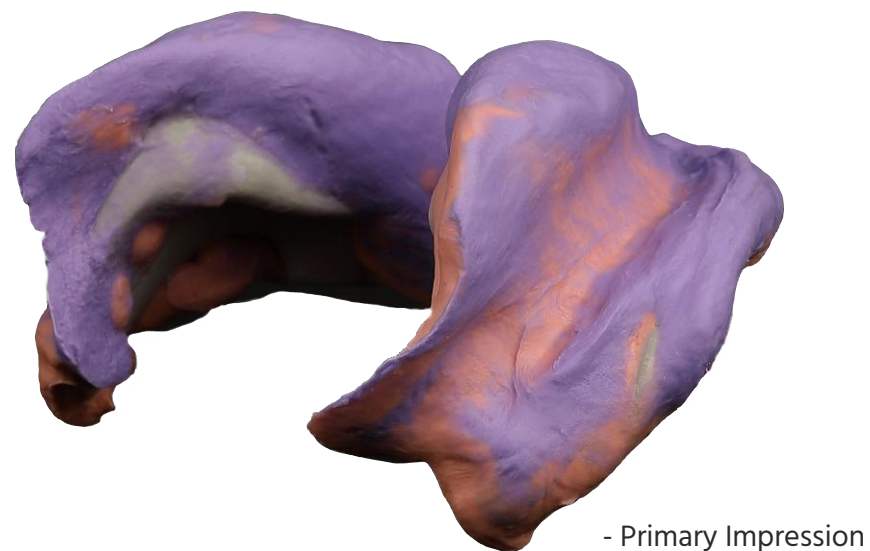

Xtend trays™ are impression trays for making primary impressions in an edentulous jaw. They are designed to follow the ridge anatomy to capture the important landmarks. Xtend trays are used with dual consistency Alginate (tray and syringe consistency) for a good impression. These trays are autoclavable at 121°C.


## Upper Trays

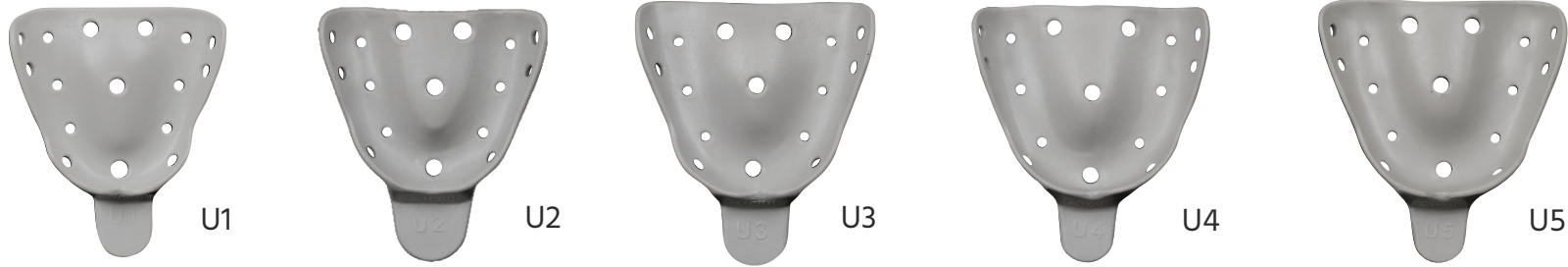


## Lower Trays - Normal Ridge



## Lower Trays - Atrophied Ridge

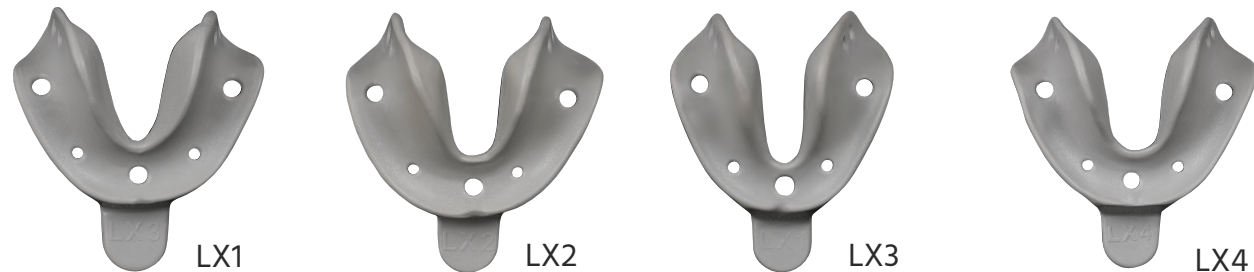


Used to select the appropriate impression trays.

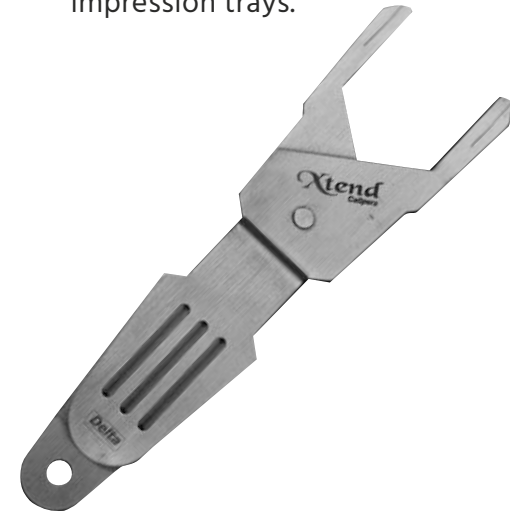

## Tray Selection:



UPPER: The distance between hamular notches are measured

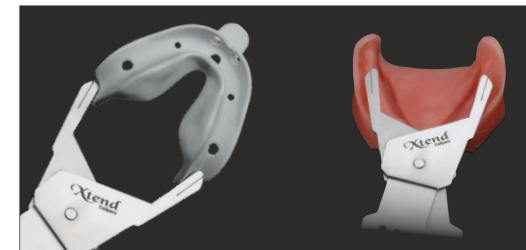


LOWER: The distance between retromolar pads are measured



Art. No. 10186

Xtend plane analyser is used to analyse the anatomical plane (Alatragal line) of the patient in relation to the interpupillary line and thereby establish the occlusal plane for the prosthesis in order to achieve functional stability.

It has a vertically adjustable (max- 10mm) interpupillary reference plate, The plane analyser is used widely in complete denture fabrication to check the occlusal plane during. Jaw Relation, Denture trial etc.

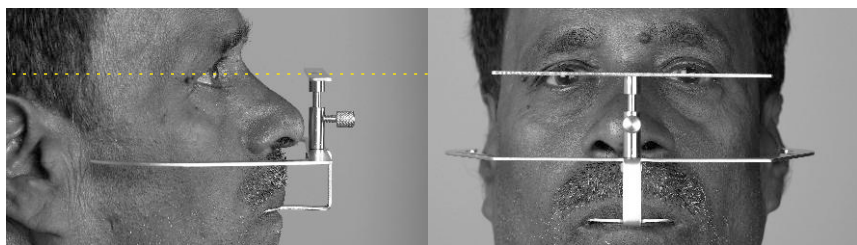



## Tentric Recorder

Tentative intervestibular relation recording tray.

Autoclavable at 121°C



Art. No. 10362

## Bio-Tracer - Intra Oral Tracer

Biotracer is a simple intra oral tracing device to determine the patient's centric and vertical dimension in an fully or partially edentulous situation.



Art. No. 10359

## Model Analysis

A hand book on detailed edentulous model analysis



Art No: 10187



## Professional Set

Art No: 10184

Content:

Upper Trays - 5 nos

Lower Trays (Normal Ridge) - 5 nos

Lower Tray (Atrophied Ridge) - 4 no

Tray Calipers - 1 no

Tentric Recorder - 1 no



## Starter Kit

Art No: 10183

Content:

Upper Trays - 3 nos

Lower Trays (Normal Ridge) - 2 nos

Lower Tray (Atrophied Ridge) - 1 no

Tray Calipers - 1 no

Tentric Recorder - 1 no



## Individual Trays

Manufactured by:

**DELTA LABS**

# 323, SIDCO Women Industrial Estate, Vellanur P.O.,

Thirumullaivoyal, Chennai - 600 062. INDIA

Ph: (+91) 9884262713 / 9884922713

Email ID: [deltalabsindia@gmail.com](mailto:deltalabsindia@gmail.com)

Website: [deltadentallabs.com](http://deltadentallabs.com)



ISO 13485:2016  
Certified

## Recommended products



Xtend Alginate Mixing Bowl

(Set of 2) - Autoclavable

Art No: 10188



Xtend Base Former

Art No:

Medium - 10180

Large - 10181



Syringe for Alginate - Autoclavable

Lascod, Italy

Art No: 801033



Tray Material : Cavex CA37\*

Art No: AA055

Syringe Material: Cavex Cream\*

Art No: AA600

\*Registered brand of Cavex Holland BV

For Training contact: +91 9941050942 [www.cadetdental.com](http://www.cadetdental.com)

